

SURGICAL MILK

Instrument Lubricant and Rust Inhibitor

- Bacteriostatic
- Clinically tested to be steam penetrable and steam sterilizable
- Effective in steam, dry or chemical vapour sterilizers
- Reduces instrument repair and replacement costs




Effective	Prevents spotting, staining and instrument rust. Keeps instruments shiny and prevents dulling.
	Will not damage your evacuation system's components. Ideal for use in both wet and dry vacuum systems.
Convenient	Not sticky, oily or greasy.
	1 part solution to 5 parts water can be used for up to 7 days.
	Can be used on a variety of metals: steel, carbide, tungsten, stainless steel and aluminum.


4 L bottle of concentrate replaces 20 L of ready-to-use solution

Use During Autoclaving on:

- Dental/Surgical Instruments
- Hinges
- Box Locks
- General Equipment

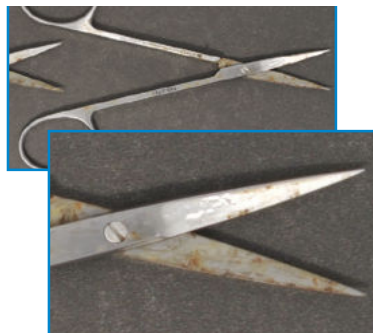

Untreated

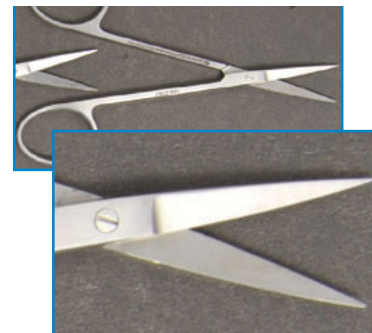

Treated with Surgical Milk



4 L
Item #:
SM-L



2 L
Item #:
SM-S



made in
Canada

To order Surgical Milk or for more information, contact your Account Manager or Germiphene at 1.800.265.9931 | info@germiphene.com | www.germiphene.com

REV 01122024