

GermiCide³

Dental & Medical, Equipment SURFACE DISINFECTANT CLEANER

- 3 Minute complete broad spectrum kill
- Cleaner & Disinfectant
- pH range 7-9
- Bactericidal, Fungicidal, Virucidal,
& Mycobactericidal (Tuberculocidal)



COVER ALL YOUR SURFACES

Dental • Medical • Janitorial • Sanitation EMS Veterinary
Pet Care Facilities Body Art Studio • Athletic Centres • Spas

SANITIZES IN 10 SECONDS

BACTERICIDAL

Staphylococcus aureus and
Enterobacter aerogenes in **10 Seconds**

- Staphylococcus aureus ATCC 6538
- Staphylococcus aureus (MRSA)
- Enterococcus faecalis (VRE)
- Pseudomonas aeruginosa ATCC 15442
- Salmonella enterica ATCC 10708
- Listeria monocytogenes ATCC 19111
- Escherichia coli ESBL
- ATCC BAA-196 (E. coli)
- Klebsiella pneumoniae (Carbapenem-resistant)
ATCC BAA-1705
- Serratia marcescens ATCC 14756 and Acinetobacter
baumannii ATCC 19606 in **2 Minutes**

FUNGICIDAL

Trichophyton interdigitale ATCC 9533 in **2 minutes**

MYCOBACTERICIDAL (TUBERCULOCIDAL)

Mycobacterium bovis BCG
(one of the causative agents of tuberculosis (TB)) in **2 Minutes**

VIRUCIDAL

Human Immunodeficiency Virus (HIV-1) in **1 minute**

- SARS-CoV-2 (USA-WA1-2020, NR-52281)
- Herpes Simplex II Virus (HSV 2) G Strain
- Influenza Virus Strain A2/Hong Kong
- Duck Hepatitis B Virus (HBV), Hepatitis C Virus (HCV)
- Avian Influenza A (H3N2) Virus
- Swine Influenza A (H1N1) Virus in 2 minutes

Rhinovirus ATCC VR-286, Polio 1 Virus ATCC VR-58
and Murine Norovirus in **3 minutes**

COMPATIBLE WITH A VARIETY OF SURFACES

- Anodized aluminum
- Brass 443
- Copper 110
- Stainless steel 304
- Stainless steel 410
- Stainless steel 420
- Acrylic
- Polystyrene
- PVC
- Neoprene
- HDPE
- Polypropylene
- Silicone
- Polycarbonate
- Vinyl

REV 01 062024


GermiCide³



- ▶ Non-flammable. No special shipping or storage is required
- ▶ Non-sticky and leaves no hazardous chemical residue when used as directed

REFILL
eco friendly
choice

Classes of Microorganisms* | GermiCide³ Disinfection

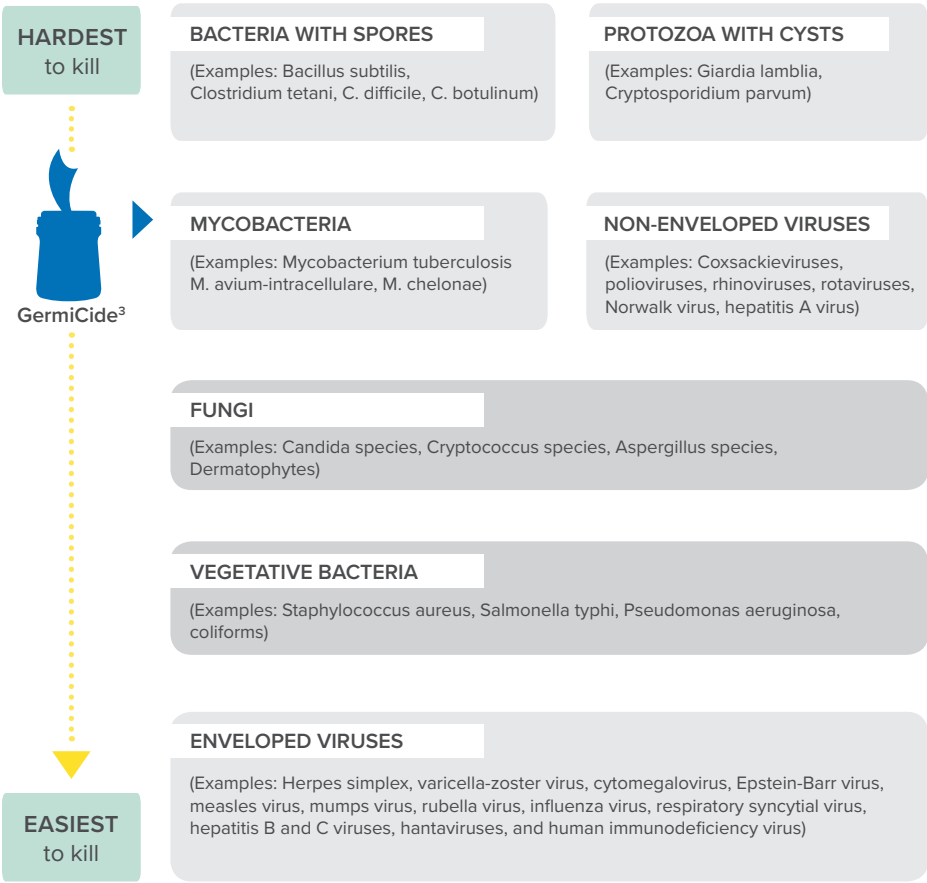

*Source: Infection Control Guidelines: Hand Washing, Cleaning, Disinfection and Sterilization in Health Care (CCDR 24S8)

Try both **liquid + wipes** to ensure all your infection control needs are met



- ACTIVE INGREDIENTS:**
- n-Alkyl dimethyl benzyl ammonium chloride.....**0.154%**
 - n-Alkyl dimethyl ethylbenzyl ammonium chloride.....**0.154%**
 - Isopropanol.....**21.000%**

To order GermiCide³ or for more information, contact your Account Manager or Germiphene at 1.800.265.9931 | info@germiphene.com | www.germiphene.com

REV 01 062024